SAMSUNG

# Mobile Excellence HM70 EVO







### **Mobile Excellence**

The HM70 EVO ultrasound system is a high-performance hand-carried ultrasound system, evolved to support a diverse range of applications. The system has streamlined workflow, durability, and high resolution imaging that can be used in a variety of clinical situations. Samsung has harmonized the key aspects of the HCU to meet today's changing needs and provide healthcare professionals the confidence they need in their work environment.





Versatility Supports various applications and



Intuitiveness

Durability



Mobility Go almost a whole working day on a single charge



Exceptional image clarity and penetration is adapted from Samsung's high-end imaging platform to produce clear, uniform, high resolution 2D and color images. The HM70 EVO delivers sophisticated imaging technologies to empower ultrasound professionals with diagnostic confidence wherever needed.

#### Noise reduction filter to improve 2D image quality

The noise reduction filter improves edge enhancement and creates sharp 2D images for optimal diagnostic performance. In addition, ClearVision provides application-specific optimization and advanced temporal resolution in live scan mode.

#### Clarify blurred area to provide clearer images

HQ-Vision<sup>™ 1</sup> provides clear images by mitigating the characteristics of ultrasound images that are slightly blurred than the actual vision.


#### Uniform imaging performance of overall image area from near-to-far

**S-Harmonic™** mitigates the signal noise, enhances contrast, and provides uniform image performance of overall image area from near-to-far.

#### **Directional power Doppler** to examine peripheral vessels

S-Flow™ uses directional power Doppler technology, enabling you to examine even the peripheral vessels. It displays information on the intensity and direction of blood flow



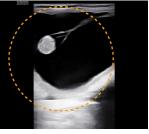

Liver with ClearVision



Shoulder




Shoulder with HQ-Vision™



Gallbladder



Gall bladder with S-Harmonic™

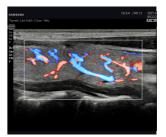

Thyroid with S-Flow™



Kidney with S-Flow™ in Dual Live

## **Increased Operational Efficiency**

The simple, intuitive interfaces enable users to immediately get to work in emergency situations. The HM70 EVO is engineered to provide the most direct interfaces and smooth actions from anywhere, avoiding delays caused by configuration. Samsung has considered even the smallest details of usability to support users with the exceptional possible performance.



#### Intuitively tuned interfaces

Intuitive buttons and screen layouts are optimized to express information clearly and enable complex actions to be simply performed.







#### View images and cines in full screen

In the Full screen mode, the ultrasound examination can be performed while viewing the image/cine that is fully expanded to the entire monitor.

#### Quickly exam even in an emergency case

In case of emergency, the **QuickExam** function allows to begin an exam with a temporary ID and veterinary information can be modified after the exam.

#### Select your transducer and preset combination with one touch

With one touch, the user can select the most common transducer and preset combinations. QuickPreset increases efficiency to make a full day of scanning simple and easy.

#### Save image data directly to USB memory

QuickSave function allows image data to be saved directly on USB memory during the exam.





Widescreen (Scan Layout 2)

#### Build predefined protocols to ensure every step is followed every time

EzExam+™<sup>1</sup> enables users to build or to use predefined protocols. It transforms the ultrasound investigation into a streamlined process. EzExam+™ ensures the full investigation is performed, eliminating the risk of forgetting an image or loop capture, as well as measurement and transducer preset changes.

#### Provide valuable reference tool to help in scanning

As a reference tool for plane scanning, **EzAssist™**<sup>1</sup>has multiple clinical applications during live scanning. This allows clinicians to practice procedures using valuable reference materials such as clinical reference images and animation.

# **Highly Usable Design**

The wave pattern of the exterior surface and the honeycomb structure of the internal design communicate strength and confidence to the healthcare professionals. The magnesium case surrounding the system and enhanced battery performance offer true usability and endurance.



#### Protect your system against damage from knocks and drops with the durable cover

Considering the importance of durability in a system that needs to be moved frequently, exterior wave pattern and interior honeycomb structure have been incorporated in the design to improve the rigidity of the casing bv 38%\*.



\* Compared to the Samsung ultrasound HM70A





#### **1** Go almost a whole working day on a single charge <sup>1,2</sup>

The battery lasts approximately 7.5 hours with the extended power pack. This is more than 100% longer than the previous model, expanding its deployment range.

#### 2 Connect up to 3 transducers and save your time from switching over transducers<sup>1</sup>

You can connect up to three transducers at once using the extended transducer ports on the optional cart, saving time and effort spend switching transducers in situations where every second counts.

#### 3 Convenient carrying case safely protects the system with shockproof cushion<sup>1</sup>

The carrying case allows you to safely and conveniently transport your whole package, including ultrasound system, transducers, and other accessories.







### Available in Various Clinical Environment

The HM70 EVO supports a wide range of clinical applications and environments — including abdomen, cardiac, vascular, musculoskeletal, emergency, and small parts. The versatility enables healthcare professionals to perform effectively on targeted examinations.



Ovary




Bladder





Wall motion measurement by Strain+<sup>1</sup>

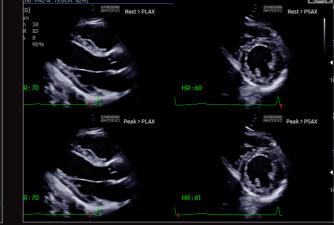

StressEcho<sup>1</sup>




Pelvis in Panoramic+<sup>1</sup>



Foal heart

		LA3-16AD/3.5cm/44Hz				_				
R 1 A 6	17 04	Mean SD	0.59 mm 0.47 mm 0.05 mm 1.00						1	
									C N	
	,	-				-		-		
		No.	-		and the second second	P HAT	100	DEP.		
				-		-	24			

IMT measurement with AutoIMT+<sup>1</sup>

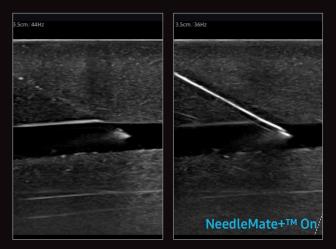

Visualization of Needle using NeedleMate+™<sup>1</sup>



Suspensory ligament

### **Comprehensive selection of transducers**

#### Curved array transducers



CA1-7AD Abdomen, Gynecology, Obstetrics, Pediatric, Musculoskeletal. Vascular. Urology, Emergency



 $C_{2-8}$ Abdomen, Obstetrics, Gynecology, Pediatric, Musculoskeletal, Vascular, Urology



CA2-8AD Abdomen Obstetrics Gynecology, Pediatric, Musculoskeletal, Vascular, Urology, Emergency

LA2-95 \*

Pediatric

Abdomen. Obstetrics.

Small Parts, Vascular,

Gynecology, Musculoskeletal,



#### CA4-10M \*

LS6-15

Pediatric

Abdomen, Obstetrics,

Small Parts, Vascular,

Gynecology, Musculoskeletal,

Pediatric, Vascular, Abdomen, Cardiac, Obstetrics, Gynecology, Musculoskeletal, Urology, Emergency, Vet Cardiac, Vet Abdomen

#### Linear array transducers



LA3-16AD Small Parts, Vascular, Musculoskeletal, Abdomen, Obstetrics, Gynecology, Pediatric, Emergency, Vet Abdomen

#### LN5-12 Small Parts, Vascular, Musculoskeletal, Abdomen. Obstetrics, Gynecology,

Pediatric

#### Phased array transducer



PA1-5A \* SP3-8 Abdomen, Cardiac, Vascular, Pediatric, TCD, Emergency

Abdomen, Cardiac, Vascular, Pediatric, TCD, Emergency, Vet Cardiac, Vet Abdomen



PN2-4 Abdomen, Cardiac, Vascular, Pediatric, TCD, Emergency



PA3-98 Abdomen, Cardiac, Vascular, Pediatric, TCD, Emergency



LA2-9R Vet \*\* Musculoskeletal. Reproduction


#### \* Ergonomic transducers

The new transducer supports natural grip by moving the max-width point to a more forward position and also increasing the length of the grip to allow balanced weight distribution.

#### \*\* Equine dedicated transducer



Cleaning and disinfection guide

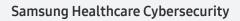

To address the emerging need for cybersecurity, Samsung provides a solution to support our customers by offering the tools to protect against cyberthreats that may compromise invaluable data and ultimately degrade the quality of care.



- \* This product, features, options, and transducers may not be commercially available in some countries. \* Sales and Shipments are effective only after the approval by the regulatory affairs.
- Please contact your local sales representative for further details. \* Some of the ultrasound images in this catalog are acquired with the human body.
- S-Vue Transducer™ is the name of Samsung's advanced transducer technology.
- \* This product is a medical device, please read the user manual carefully before use.

1. Optional feature which may require additional purchase.

2. It is assumed that the average working hours per day is 8 hours. This is based on internal study and compared to the previous model, HM70A The 7.5 hours of battery life is the combined battery life of the system and extended power pack

### SAMSUNG MEDISON CO., LTD.

© 2023 Samsung Medison All Rights Reserved. Samsung Medison reserves the right to modify the design, packaging, specifications, and features shown herein, without prior notice or obligation.